

# 66 SERIES FFC

## High Flow Liquid and Gas Housing



The **66 Series FFC** housing is a version of Amazon's successful **66 Series** range and specifically designed for use with the large bore **SupaPleat FFC** and **SupaSpun FFC** cartridges. This provides a filtration system that offers maximum efficiency, high flow rates and long life. The high quality electropolished 316L stainless steel housing is easy to use and is suited to many filtration processes across a range of industries.

They are available in a selection of sizes and with a choice of connections and closure types, therefore the optimum housing can be provided for each application.





The **66 Series FFC** housing is available to be used with 20", 40" or 60" cartridges. The cartridge is normally operated in 'out to in' mode.

Housings can be installed in either the vertical or horizontal position and can be supplied with fixed or adjustable legs. This ensures that the housing can normally be used even where space or changes to inlet and outlet pipework are limited.

### Product Features

- Designed for high flow systems
- Easy and quick cartridge change-out
- Electropolished inside and out
- Available for 20", 40" or 60" filters
- Each housing is hydrostatically pressure tested
- Available with 3" or 4" connections and in a range of types
- Can be supplied with adjustable or fixed legs where required
- Rated at 10 bar g liquid, 3 bar g gas
- PED and ATEX Compliant

### Accessories

Pressure Gauge and Vent Valve Kit	
Pressure Gauge and Auto Vent Valve Kit	
Ball Valve and Nipple	
Cast Vee Clamp with Interlock (3 and 5 stack only)	



## Features and Benefits

- Designed for use with the **SupaPleat FFC** and **SupaSpun FFC** cartridges for maximum flow rates
- Utilises large bore cartridge allowing quick and easy filter change out
- Can be used vertically or horizontally with or without legs
- Flexible design allows use of either V band or Rathmann clamp closure
- Full housing customisation service available
- Unique safety interlock device available to prevent accidental opening under pressure

## Product Range

Amazon Filters design and manufacture an extensive range of liquid filtration systems, providing one of the widest selections of filtration vessels, depth filters and pleated cartridges available today, ensuring excellent process and filter compatibility. As one of Europe's leading manufacturers of filtration systems, Amazon Filters is able to offer a solution to meet all filtration requirements. Amazons' filters are manufactured in clean room conditions and all products come with detailed Product Validation Guides and Technical Support documentation.



**AMAZON**

# 66 Series FFC Technical Data

## Maximum Operating Conditions

Pressure Range:	Liquid:	Vee clamp:	10 bar g
		Bolted:	FV / 10 bar g
	Gas:	Vee clamp:	3 bar g
		Bolted:	FV / 3 bar g
Temperature Range:	-10°C / 150°C		

## Materials of Construction

Body:	316L St.Steel
Closure:	Vee Clamp - 304 St. Steel
	Bolts - Zinc plated carbon steel
	(Optional stainless steel)

## Finish

Vessel:	Electropolished
---------	-----------------

## Connections

Inlet/Outlet:	3" Flanged
Vent:	Horizontal: ½" BSPM
	Vertical: ¼" BSP
Drains:	½" BSPM

Model	Dimensions					Volume (Litres)	Dry Weight (kg)
	A	B	C	D	E		
66012HAC30C1EE	1070	853	214	201	550	27.4	31.0
66014HAC30C1EE	1599	1382	214	201	1000	46.6	36.3
66014HAI30C1EE	1599	1382	214	201	1000	46.6	40.0
66016HAI30C1EE	2102	1885	214	201	1550	59.1	45.0

## PED Status:

Fluid Groups:	Compliant with 2014/68/EU
	Liquids Group 1 + 2
	Gas Group 2


## Marking:

Exceptions:	CE
	Group 1 liquids - Explosive

## ATEX:

Compliant with 2014/34/EU, level CE Ex II 2 GD c X

Technical drawings for each individual housing showing dimensions, volume, weights, maximum temperatures and design pressures available on request.

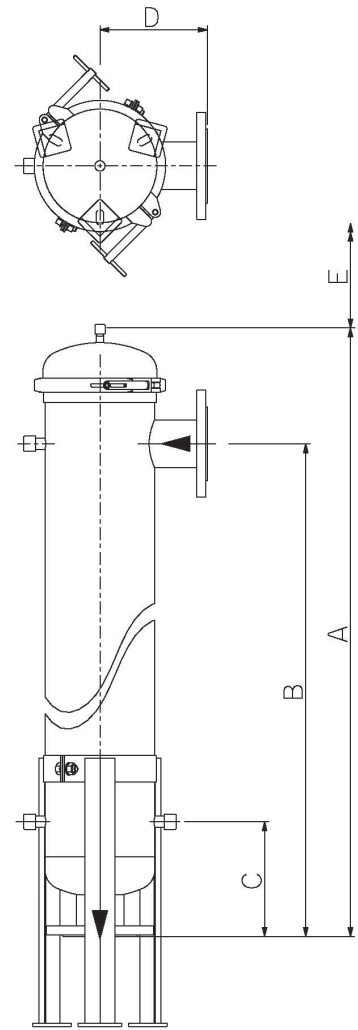

# Ordering Guide

6601	4	H	A	I	30	C	1	E	E
No of Stacks	Cartridge Length	Cartridge End Cap Style	Closure Style	Mounting	Connection	Connection Type	Materials	Finish	Seal
6601 - 1 (Available up to 30 round - please contact the Sales Office)	2 - 20" 4 - 40" 6 - 60"	H - FFC Cartridge	A - Cast Vee Clamp H - Hinged Plus Bolts U - Cast Vee Clamp with Interlock W - Hinged + St.St Bolts	I - Adjustable Legs - Horizontal (40" and 60" only) C - Adjustable Legs - Vertical (20" and 40" only)	30 - 3" 40 - 4"	A - BSPM C - EN1092-1 Flange PN16 G - Table 'E' Flange H - ASME Flange 150# K - Stub Pipe L - BS4825-3 Tri Clamp M - DIN11851 N - RJTM P - SMS M (Swedish) R - IDFM S - DIN2633 Flange	1 - 316 st.st (L Grade)	E - Electropolished	E - EPDM F - FEP / Silicone (Flat cover only) S - Silicone V - Viton

Example: 66014HAI30C1EE= 66 Series, 1 x 40" FFC cartridge, Cast vee clamp closure, Adjustable legs - Horizontal, 3" EN1092-1 Flange PN16 connections. 316L St.Steel construction with EPDM gaskets, electropolished all over.

# AMAZON FILTERS LTD.

Albany Park Estate, Camberley, Surrey, GU16 7PG, ENGLAND

Tel: +44 (0) 1276 670 600 Email: sales@amazonfilters.co.uk Web: www.amazonfilters.com

66FFC - Issue 05